


## ABSTRAK

Ada dua sistem dalam dental unit, yaitu sistem hidrolik dan sistem motorulir. Untuk sistem hidrolik, setiap tidak dipakai, kursi diturunkan pada posisi terendah, dan posisi backrest rebah maksimal.

Ini dilakukan karena ketika di posisi tinggi, piston terbuka sehingga memungkinkan kotoran menempel di piston, dan ketika piston dipakai, kotoran akan masuk ke sela-sela piston, menyebabkan o-ring menjadi aus. Jika terjadi kebocoran, cek ketersediaan oli, kalau oli habis, silakan diisi lagi. Sebaiknya piston diberi pelumas setiap sebulan sekali.

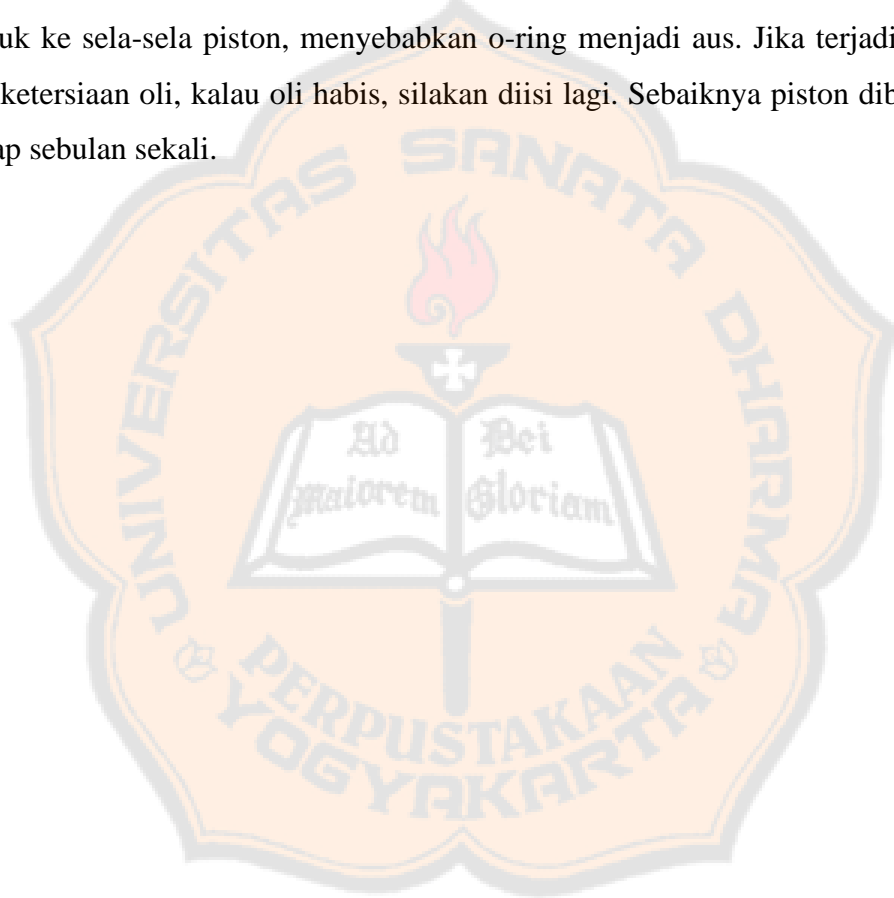

## ABSTRACT

*There are two systems in the dental unit, namely the hydraulic system and the screw motor system. For the hydraulic system, every time it is not used, the seat is lowered to its lowest position, and the maximum backrest position.*

*This is done because when it is in the high position, the piston is open allowing dirt to stick to the piston, and when the piston is worn, dirt will enter between the pistons, causing the o-rings to wear out. If there is a leak, check the availability of oil, if the oil runs out, please refill it. It is better if the piston is lubricated once a month.*

